




# ISO CERTIFICATE



# CE CERTIFICATE



ADD: Dingtao industrial area shandong province china.  
 PC:274199  
 Tel:86-0530-2135555  
 fax:86-0530-2135559  
 e-mail:sdwuzhouyx@163.com  
 Web:https:infusionsetfactory.com  
<https://sdwuzhoumedical.en.made-in-china.com>

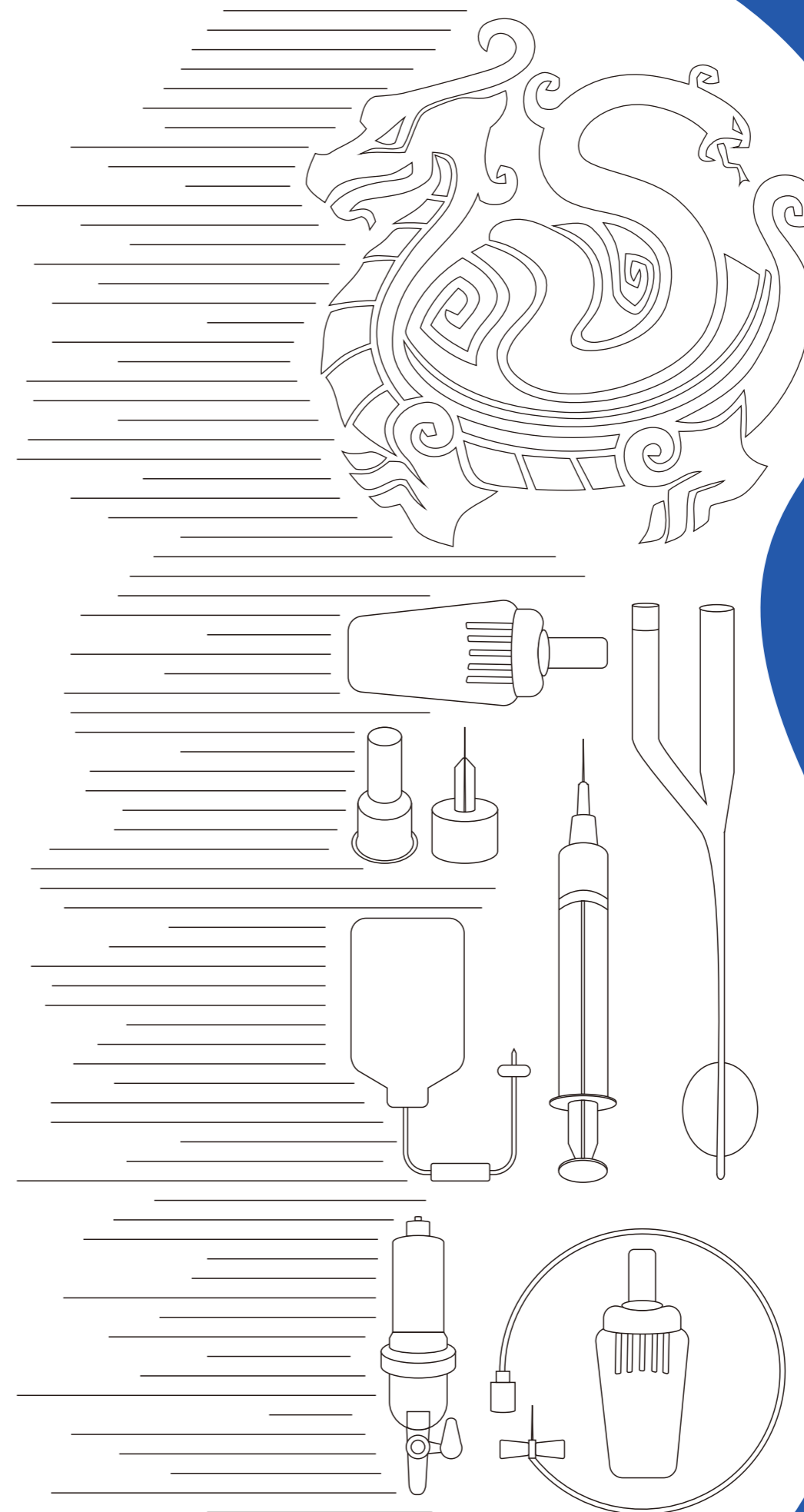

# 产品画册 PRODUCT CATALOGUE

山东五洲医疗器械有限公司  
 SHANGDONG WUZHOU MEDICAL  
 EQUIPMENT CO.LTD





山东五洲医疗器械有限公司  
SHANGDONG WUZHOU MEDICAL  
EQUIPMENT CO.LTD



天津市麦迪克医疗器材有限公司  
TIANJIN MEDIC MEDICAL  
EQUIPMENT CO.LTD

## 企业简介

Shandong Wuzhou Medical Equipment Co., Ltd, is strategically located in Heze city, Shandong Province, with proximity to the Qingdao Port. Over the past decades, our company has dedicated itself to the production of high-end disposable medical equipment products. Committed to precision in market positioning, technical innovation, and seamless product quality, our company has a strong technical foundation, dedicated research and development, and a skilled production team with access to top-notch quality testing instruments.

### Main products:

Syringe,	Insulin Syringe,
Infusion Set,	Insulin Pen Needle,
Blood Transfusion Set,	Peripheral Blood,
Burette Type Infusion Set,	Dental Needle,
Feeding Syringe,	Latex Urinary Catheters.

We consistently adhere to ISO 13485 standards and European CE directives, holding registrations in numerous countries. Our products have achieved widespread market penetration, reaching customers across Asia, Africa, Europe, Latin America, and more than 30 other countries and regions. Our unwavering focus on product quality and excellent after-sales service has garnered praise from both domestic and international customers.



# + CONTENTS

**01** 采血针系列  
BLOOD TAKING NEEDLE

**02** 导尿管系列  
CATHETER

**03** 吊瓶系列  
BURETTE SET

**04** 输血器系列  
BLOOD TRANSFUSION SET

**05** 输液器系列  
INFUSION SET

**06** 胰岛素系列  
INSULIN

**07** 头皮针系列  
SCALP NEEDLE SERIES

**08** 注射器系列  
SYRINGE SERIES

**09** 牙科针系列  
DENTAL NEEDLE

**10** 注射针系列  
INJECTION NEEDLE

**11** 配件系列  
ACCESSORIES

采血针系列  
BLOOD TAKING NEEDLE



产品特点 Product Features

- 1、集采血笔和刺血器于一体的双重结构；
- 2、安全、避免交叉感染；
- 3、符合人体工程学结构，操作简单；
- 4、独特的结构设计，避免二次使用。

- 1. Dualstructure integrating blood collection pen and lancet;
- 2. Safe and avoiding cross-infection;
- 3. Ergonomic structure and simple operation;
- 4. Unique structural design to avoid secondary use.

产品参数 Product Parameters

Disposable Blood Taking Lancet			
Specs	Min. Outer Diameter	Max. Outer Diameter	Pantone Color
0.2 (33G)	0.203	0.216	113C
0.23 (32G)	0.229	0.241	3275C
0.25 (31G)	0.254	0.267	306C
0.3 (30G)	0.298	0.320	Puple C
0.36 (28G)	0.349	0.370	286C
0.45 (26G)	0.440	0.470	355C
0.6 (23G)	0.600	0.673	1235C

导尿管系列  
CATHETER

Disposable Sterile Catheter  
一次性使用无菌导尿管

Made of natural latex material, mainly composed of excretion cone interface, tube body, eyelet, tip, balloon (singlecavity without capsule, double cavity and above with capsule), valve, filling cone interface, flushing cone Theinterface is composed of products that are sterile after ethylene oxide sterilization.

由天然乳胶材料制成，主要由排泄锥形接口、管身、孔眼、尖部、球囊（单腔为无囊，双腔及以上有囊）、阀、充起锥形接口、冲洗锥形接口组成，经环氧乙烷灭菌产品无菌。



产品参数 Product Parameters

Model Specs	球囊容积				
	single lumen	double lumen infant use	double lumen female type	double lumen/men type(standard)	three-cavity type
2.0mm(6Fr/Ch)	/	3ml, 3~5ml, 5ml	/	/	/
2.7mm(8Fr/Ch)	/	3ml, 3~5ml, 5ml	/	/	/
3.3mm(10Fr/Ch)	/	3ml, 3~5ml, 5ml	/	/	/
4.0mm(12Fr/Ch)	/	/	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	/
4.7mm(14Fr/Ch)	/	/	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	/
5.3mm(16Fr/Ch)	/	/	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml, 30ml, 40ml, 30~50ml, 50ml, 60ml, 80ml
6.0mm(18Fr/Ch)	/	/	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml, 30ml, 40ml, 30~50ml, 50ml, 60ml, 80ml
6.7mm(20Fr/Ch)	/	/	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml, 30ml, 40ml, 30~50ml, 50ml, 60ml, 80ml
7.3mm(22Fr/Ch)	/	/	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml, 30ml, 40ml, 30~50ml, 50ml, 60ml, 80ml
8.0mm(24Fr/Ch)	/	/	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml, 30ml, 40ml, 30~50ml, 50ml, 60ml, 80ml
8.7mm(26Fr/Ch)	/	/	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml, 30ml, 40ml, 30~50ml, 50ml, 60ml, 80ml
9.3mm(28Fr/Ch)	/	/	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml, 30ml, 40ml, 30~50ml, 50ml, 60ml, 80ml
10.0mm(30Fr/Ch)	/	/	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml 30ml, 40ml, 30~50ml, 50ml	10ml, 5~10ml, 15ml, 5~15ml, 20ml, 30ml, 40ml, 30~50ml, 50ml, 60ml, 80ml



吊瓶系列  
HANGING BOTTLE



一次性使用滴定管式输液器

Disposable Burette Infusion Set

- (1) 带切断阀的100ml和150ml刻度滴定管
- (2) 60滴=1±0.1ml
- (3) 带过滤器的柔性滴水室
- (4) 用保护帽滑动或锁定适配器
- (5) 乳胶/无乳胶注射部位
- (6) ABS钉，带盖通风，带过滤器的滴水室，带Y注入点，鲁尔锁接头，管道1.5m/1.8m/2m
- (7) 医用级PVC管，可提供非DEHP PVC
- (8) 医用TPEG滴定管

可根据客户要求选择:

1. 带或不带排气塞;
2. 带或不带Y注入端口
3. 鲁尔卡瓦或鲁尔锁
4. 有无后眼
5. 21-25G针

- (1) 100ml and 150ml graduated burette with shut-off valve
- (2) 60 drops=1±0.1ml
- (3) Flexible drip chamber with filter
- (4) Slip or lock adapter with protective cap
- (5) Latex/Latex-free injection site
- (6) ABS Spike, air vented with cap, drip chamber with filter, with Y injection site, luer lock connector, tubing 1.5m/1.8m/2m
- (7) Medical grade PVC tube, NON-DEHP PVC available
- (8) Burette made by medical TPEG

Optional for customer requirement:

1. With or without air vented spike;
2. With or without Y injection port
3. Luer slip or luer lock
4. With or without backeye
5. 21-25G needle

输血器系列  
BLOOD TRANSFUSION SET



一次性输血器

Disposable Transfusion Set

包装细节:

- (1) 带或不带通风口的钉子
  - (2) 8.8厘米滴水室，配备62毫米17.3微米过滤器
  - (3) ABS轭夹
  - (4) 乳胶球
  - (5) 鲁尔滑动连接器
  - (6) 鲁尔帽或带针
- 用于血液或血液成分的输血。  
模制乳胶冲洗球注射膜，用于间歇用药。  
大约20滴/毫升。
- (7) 由医用级PVC制成的透明软管。
  - (8) 带200微米过滤器和18克皮下注射针

包装:PE包装或纸塑包装。

可用选项

- 有或没有通风口。有没有针
- 有或没有“Y”部位
- 乳胶或无乳胶
- 鲁尔锁或鲁尔滑

尺寸:

18G/19G/20G/21G/22G/23G/24G/25G/26G/27G

Packaging Details :

- (1) Spike with or without air vent.
  - (2) 8.8cm drip chamber with 62mm 17.3micron filter
  - (3) ABS roller clamp
  - (4) Latex bulb
  - (5) Luer slip connector
  - (6) Luer cap or with needle
- Designed for transfusion of blood or blood component.  
Moulded latex flush ball injection membrane for intermittent medication.  
Approximately 20 drops/ml.  
Clear, transparent and soft tubing manufactured from
- (7) Medical grade PVC.
  - (8) With 200 micron filter and 18 G hypodermic needle.


Packing: PE package or paper-plastic package.

Options available.

- With or without Airvent.
- With or without Needle
- With or without “Y” site(Latex or Latex-free)
- Luer lock or luer slip

Size :

18G/19G/20G/21G/22G/23G/24G/25G/26G/27G

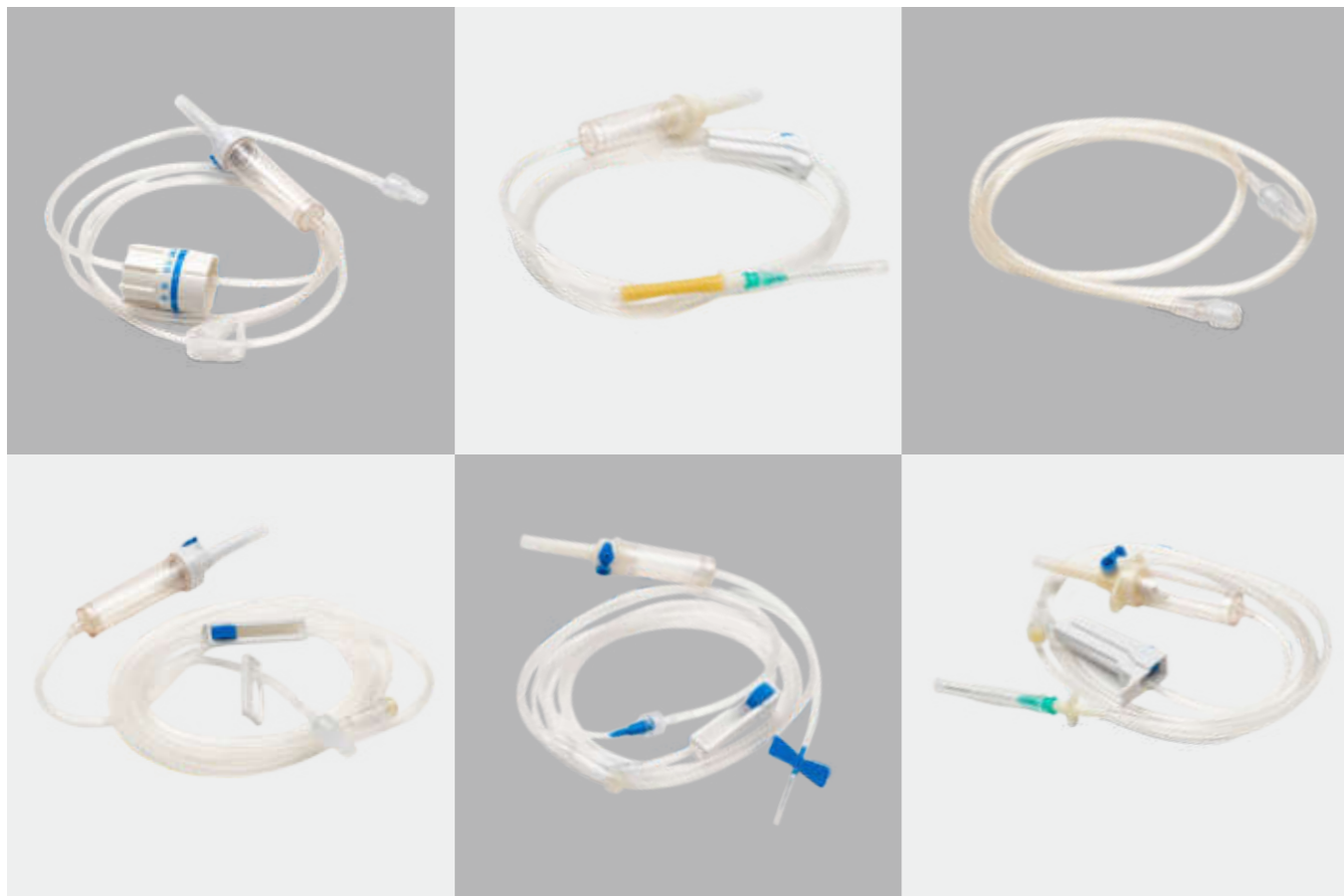

产品参数 Product Parameters

Accessories	Description
Commodity	Disposable iv infusion set
Material	PVC tube, PE chamber
Sterile	Eo gas sterile, non-toxic, non-pyrogenic
Type	10 drops/ml or 15 drops/ml or 20 drops/ml or 60 drops/ml
Tube	120cm, 135cm, 150cm, 180cm, 200cm, 250cm
Certificate	CE&ISO13485
MOQ	20000PCS
Supply Ability	700000PCS per day
Scope of Application	Used in the vein drop for the human body
Attention	Discard after use; For single use only



滴液室：分为48mm、54mm、60mm、62mm滴斗  
通风钉：带或不带通风口、鲁尔锁或直口  
流量调节器：由聚乙烯或ABS制成  
柔性软管：1.5m软PVC，抗扭结  
保护帽：用于聚乙烯制成的带内螺纹的封闭穿孔装置，防止细菌进入  
“Y”注入端口：带或不带  
橡胶注射部位：有无

**Drip Chamber** : the drip chamber can be divided into 48mm,54mm,60mm,62mm.  
**Vented Spike** : with or without air vent Luer lock or luer slip  
**Flow Regulator** : made of polyethylene or ABS  
**Flexible Soft Tube** : 1.5m soft PVC and kink resistant  
**Protective Cap** : for closure piercing device made of polyethylene with internal thread that prevents the bacteria from coming in  
**“Y” Injection Port** : with or without  
**Rubber injection site** : with or without

胰岛素系列  
INSULIN



一次性使用胰岛素注射器

Disposable Insulin Syringe

无乳胶、无毒、无热原  
几乎无痛且舒适的注射  
粗体、清晰、数字大，便于阅读  
E.O气体灭菌  
5年保质期0.5m  
规格：0.3ml 0.5 ml 1 ml



Disposable Insulin Syringe

Latex free, Non-toxic, Non-pyrogenic  
Nearly painless injection with comfort  
Bold, clear and large number for easy readability  
E.O Gas sterilization  
5-year shelf life  
Specification:0.3ml 0.5 ml 1 ml

产品参数 Product Parameters

Properties	Injection&Puncture-Instrument	Place of Origin	QingDao China
Brand name	Medtpoint,CAREMED	Name	Insulin syringe 1ml,0.5ml,0.3ml
Instrument classified	Class	Needle	29G 30G 31G 32G 33G 34G
Color	U-100 U-50 U-30	Fearture	Thin-walled technoloyp
Certificate	CE/ISO/FDA	Sterilization	EO gas
Customized Service	OEM/ODM/OBM	Shelflife	5 years
Sample	Free sample	Package	Individual pe bag/blisterpe 100pcs/box,20boxes/ctn

头皮针系列  
SCALP NEEDLE SERIES

单叶片 Scalp Vein Set-Single Wing

**产品结构：**由针座、针管、保护套组成。

**产品用途：**适用于对人体皮内、皮下、肌肉、静脉等注射

**产品特点：**  
针管采用SUS304不锈钢制造；  
针管薄壁设计、内径大、流量高；  
针尖三斜面结构，锋利度好；  
圆锥接头按6:100标准设计，与医疗器械适配；  
安全壳上设有防滑纹路，方便单手操作，有效防止针刺伤  
避免血液溢出感染，保障临床人员的安全。

**Product structure :** it is composed of needle base, needle tube and protective sleeve.

**Application :** suitable for intradermal, subcutaneous, intramuscular and intravenous injection of human body.

**Product features :**  
the needle tube is made of SUS304 stainless steel;  
thin wall design of needle tube, large inner diameter and high flow;



双叶片 Scalp Veinset-Butterfly

Metric	British Measurement	Color Code	Regular Specs	Packing Quantity
0.40	27G	medium gray	0.40*13mm	2000pcs/box
0.45	26G	brown	0.45*16mm	2000pcs/box
0.50	25G	ornage	0.50*16mm	2000pcs/box
0.55	24G	medium purple	0.55*20mm	2000pcs/box
0.60	23G	dark blue	0.60*25mm	2000pcs/box
0.70	22G	black	0.70*32mm	2000pcs/box
0.80	21G	dark green	0.80*38mm	2000pcs/box
0.90	20G	yellow	0.90*38mm	2000pcs/box
1.10	19G	creamy	1.10*38mm	2000pcs/box
1.20	18G	pink	1.20*38mm	2000pcs/box



注射器系列  
SYRINGE SERIES



一次性注射器  
Disposable Syringe

符合国家标准GB15810-1995和国际标准ISO7886-1。  
类别：锁定式、直进式  
形式：中头型和部分头型  
规格：1ml, 2ml, 5ml, 10ml, 20ml, 30ml, 50ml。

It conforms to national standard GB15810-1995 and international standard ISO7886-1.  
Category : luer lock , luer slip  
Form : middle-head type and partial-head type  
Specification :  
1ml, 2ml, 5ml, 10ml, 20ml, 30ml, 50ml.

产品概述

本产品主要以医用级PVC为原料，经过注塑、吹塑、组装、包装、灭菌等工序，在10万级净化车间生产。本产品采用环氧乙烷灭菌，无菌，无热原，安全可靠，使用方便。

Product Overview

This product is mainly made of medical-grade PVC as raw material, and is produced in a 100,000-level purification workshop through injection molding, blow molding assembly, packaging, sterilization and other processes. This product is sterilized by ethylene oxide, sterile, non-pyrogenic, safe, reliable and easy to use.



牙科针系列  
DENTAL NEEDLE



SHANGDONG  
WUZHOU MEDICAL


注射针系列  
INJECTION NEEDLE

注射针用的针管由我们公司自行生产;  
精准的模具设计和先进的注塑工艺;  
产品组装由先进的全自动注射针组装机完成;  
精密的生产工艺使产品刺入时的疼痛非常小。

The needle tube for injection needle is produced by our company;  
Accurate mold design and advanced injection molding technology;  
The product assembly is completed by the advanced automatic injection needle  
assembly machine:  
The precise production technology makes the pain of product penetration very small.



针号示意图  
NEEDLE NUMBER DIAGRAM

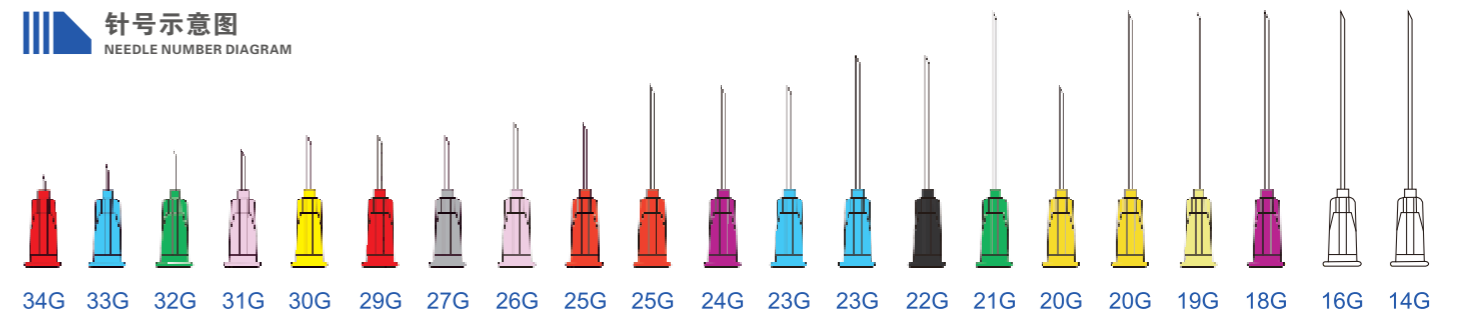



配件系列  
ACCESSORIES



公司实景  
Company Reality

